

# MEDLANCE<sup>®</sup> *plus*

## Safety Lancets

### Product features

Comfort  
MAX

CE 0344

- **Safe:** single use, self-destructing. No possibility of reuse.
- **Simple:** intuitive 2 step operation. Blood sampling simpler than ever before.
- **Sterile:** needles are sterilized with Gamma irradiation.



- **Gentle:** ultra sharp needles, ideally positioned during skin penetration.
- **Tailored:** available in 5 different versions to meet specific blood sample needs and address the variety of patients' skin types.
- **Efficient:** a "technique-independent" design provides an adequate blood flow.



## BOX

100ct 200ct 50ct 25ct 12ct

# MEDLANCE<sup>®</sup> *plus*

## Medlance<sup>®</sup> *plus SuperLite*

**30G**  
Needle Tip

1.2 mm

Cat. No. 7241, Sales Box 200ct  
Cat. No. 7242, Sales Box 100ct  
Cat. No. 7286, Sales Box 50ct  
Cat. No. 7276, Sales Box 25ct  
Cat. No. 7053, Sales Box 12ct



**Comfort**  
**MAX**

Collecting blood has never been so simple! Specially designed 30G needle tip provides higher patient's comfort without compromising usability.

## Medlance<sup>®</sup> *plus Lite*

**25G**  
Needle

1.5 mm

Cat. No. 7043, Sales Box 200ct  
Cat. No. 7243, Sales Box 100ct  
Cat. No. 7287, Sales Box 50ct  
Cat. No. 7277, Sales Box 25ct  
Cat. No. 7057, Sales Box 12ct



Ideal for blood glucose monitoring purposes.

## Medlance<sup>®</sup> *plus Universal*

**21G**  
Needle

Penetration Depth  
1.8 mm

Cat. No. 7044, Sales Box 200ct  
Cat. No. 7244, Sales Box 100ct  
Cat. No. 7288, Sales Box 50ct  
Cat. No. 7278, Sales Box 25ct  
Cat. No. 7058, Sales Box 12ct



Whenever you need larger blood samples for testing glucose, hemoglobin, cholesterol, blood type studies, coagulation, blood gas and many others.

## Medlance<sup>®</sup> *plus Extra*

**21G**  
Needle

Penetration Depth  
2.4 mm

Cat. No. 7045, Sales Box 200ct  
Cat. No. 7245, Sales Box 100ct  
Cat. No. 7289, Sales Box 50ct  
Cat. No. 7279, Sales Box 25ct  
Cat. No. 7059, Sales Box 12ct



Patient's skin is tough? More blood is needed? Reach for Medlance<sup>®</sup> plus Extra and collect larger blood sample also from men and physical workers.

## Medlance<sup>®</sup> *plus Special*

**0.8mm**  
Blade

Penetration Depth  
2 mm


Cat. No. 7046, Sales Box 200ct  
Cat. No. 7246, Sales Box 100ct  
Cat. No. 7285, Sales Box 50ct  
Cat. No. 7280, Sales Box 25ct  
Cat. No. 7771, Sales Box 12ct



Medlance<sup>®</sup> plus Special is perfect for tests requiring larger blood samples. Ultra-sharp blade provides adequate blood volume.

## Instructions for use

### 1. Twist

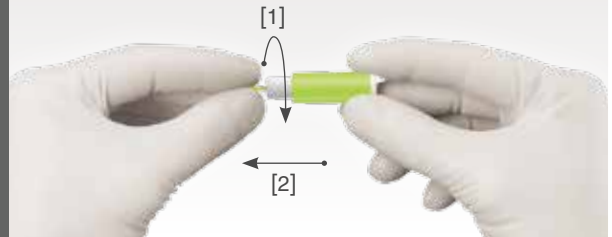

Twist off protective cap and then pull it straight out.

### 2. Press



Press the Medlance<sup>®</sup> plus body firmly against the puncture site to activate the device. Dispose of used lancets in an approved sharps container.

### 3. Collect



Gently apply intermittent pressure near the puncture site to obtain the required blood volume.

**HTL** STREFA  
high tech lab

**Distributed in the US by:**  
**HTL-STREFA, Inc.**

3005 Chastain Meadows Pkwy, Suite 300,  
Marietta, GA 30066, USA  
T: +1 770 528 0410, F: +1 770 528 0411

[info@htl-strefa.com](mailto:info@htl-strefa.com) • [www.htl-strefa.com](http://www.htl-strefa.com)

**Manufactured by:**  
**HTL-STREFA S.A.**

Adamówek 7,  
95-035 Ozorków, Poland  
T: +48 22 587 05 72, F: +48 22 811 05 36  
[info@htl-strefa.pl](mailto:info@htl-strefa.pl) • [www.htl-strefa.pl](http://www.htl-strefa.pl)